

ENDOSAFE® - PTS CARTRIDGES

SINGLE-TEST, DISPOSABLE CARTRIDGES FOR ENDOTOXIN DETECTION

FOR USE WITH ENDOSAFE® PORTABLE TEST SYSTEM AND MULTI-CARTRIDGE SYSTEM READERS

INTENDED USE

Disposable test cartridges contain Limulus Amebocyte Lysate (LAL), endotoxin, and synthetic color-producing substrate. They are intended for use with Endosafe® Portable Test System (PTS) to perform quantitative detection of endotoxins by kinetic chromogenic methods.

CAUTION:

This product is not licensed by the FDA. The Endosafe® - PTS cartridges are not to be used for end product release of Pharmaceutical Drugs, Devices or Biologics that are regulated by the FDA.


BACKGROUND AND SUMMARY

Frederick Bang observed that bacteria causes intravascular coagulation in the American Horseshoe Crab, *Limulus polyphemus*¹. In collaboration, Levin and Bang⁵ found that the agent responsible for the clotting phenomena resides in the crab's amebocytes, or circulating blood cells, and that pyrogens (bacterial endotoxin) triggered the enzymes involved in the clotting cascade.

The LAL test is the most sensitive and specific means available to detect and measure bacterial endotoxin, a fever-producing byproduct of gram-negative bacteria, commonly known as pyrogen. The basis of the test is that the endotoxin produces changes in the appearance of LAL that are easily measured.^{1-4,5} The simplicity and overall economy of the LAL Test encourages the testing of in-process solutions and raw materials as well as end-product drugs, devices, and biologics¹⁰. The USP Bacterial Endotoxins Test and the USFDA Guideline for LAL testing provide standard methods for validating the LAL Test as a replacement for the rabbit pyrogen test.^{7,9}

BIOLOGICAL PRINCIPLES


In this assay, bacterial endotoxins initiate the activation of a cascading series of serine proteases in LAL. The last activated enzyme in this series, the pro-clotting enzyme, cleaves a peptide from an endogenous substrate called coagulogen. The modified substrate produces an opacity and gelatin in LAL that is easily detected. A synthetic analog to coagulogen can also be used to quantitatively measure the endotoxin mediated activation of the LAL pro-clotting enzyme. This synthetic substrate undergoes cleavage, resulting in the release of the chromophore, p-nitroaniline (pNA). PNA is a yellow color that is measured photometrically at 385-410 nm. With the aid of a spectrophotometer, a kinetic colorimetric assay may be done, in which the early onset of color can be detected and precisely measured.

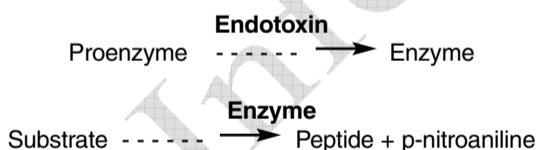


The Endosafe® - PTS cartridge and its interface with the reader have been designed to mimic currently licensed Endpoint Chromogenic and Kinetic Chromogenic Methods by measuring color intensity directly related to the endotoxin concentration in a sample. Each cartridge contains precise amounts of FDA licensed LAL formulations, chromogenic substrate, and Control Standard Endotoxin (CSE).

CARTRIDGE REAGENTS


Each Endosafe® - PTS cartridge contains four channels to which LAL reagent and a chromogenic substrate have been applied. Two of the four channels also contain an endotoxin spike, and serve as the positive product controls. The use of two Sample Channels and two Spike Channels is consistent with current FDA licensed quantitative Limulus Amebocyte Lysate (LAL) methods. The PTS reader can be used to detect endotoxin in variety of ranges.

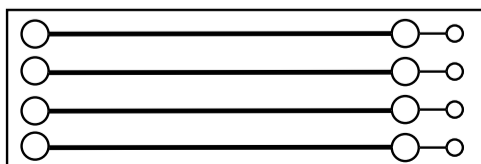

4 Channel PTS Cartridge

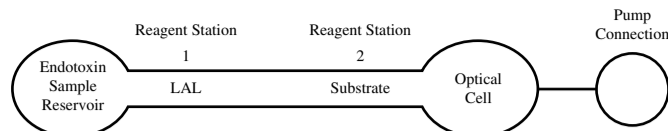

Sample Channel Zoom

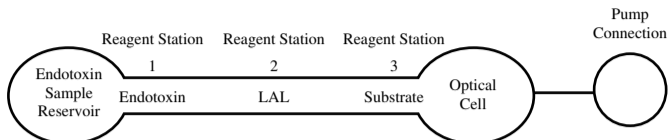

Spike Channel Zoom

STORAGE CONDITIONS AND PRECAUTIONS

PTS cartridges are relatively heat stable and should be stored at 2-25 °C. Allow the cartridges to come to room temperature before opening the pouch and testing. Prolonged exposure to temperatures above 25 °C should be avoided. To minimize contamination of the sample wells, the cartridge should be used immediately once the foil pouched seal has been opened. Cartridges are for single-test use only.

REAGENTS REQUIRED BUT NOT SUPPLIED

LAL Reagent Grade Water must be used during the initial qualification of each lot of PTS cartridges.

MATERIALS REQUIRED BUT NOT SUPPLIED

Pipettor (Endosafe® PTS 400 or equivalent) and sterile tips.

Disposable, endotoxin-free glass dilution tubes or sterile, disposable polystyrene tubes (Endosafe® T300 or equivalent) for sample collection or dilution if necessary.

Vortex-Type Mixer (if necessary).

EQUIPMENT REQUIRED BUT NOT SUPPLIED

Endosafe® Portable Test System (PTS) Reader: The reader is a dedicated instrument that accepts the cartridge and runs the PTS LAL test. The reader consists of an incubating chamber, a sample pump, four LEDs and four detectors, an alphanumeric key pad with built-in LCD, and a microprocessor. The reader operates using standard AC power or an internal rechargeable battery. Battery power also acts as automatic backup power in case of AC power failure.

Endosafe® Multi-Cartridge System (MCS): This reader is a dedicated instrument that is equipped with five independent cartridge readers. Each MCS instrument can test up to five cartridges. The cartridges can be run simultaneously or in a random access fashion. Each MCS cartridge reader consists of an incubating chamber, a sample pump, four LEDs, and four detectors.

SPECIMEN COLLECTION AND PREPARATION

Specimen for testing must be collected and prepared using depyrogenated materials and endotoxin-free reagents. Glassware must be depyrogenated by validated conditions, such as 30 minutes exposure at 250°C.⁹ It is prudent to test for endotoxin those materials that cannot be heat sterilized or those which are sold without an endotoxin-free label. Use aseptic technique at all times.

PERFORMANCE CHARACTERISTICS

Acceptance criteria for a valid assay is an archived curve correlation coefficient of ≤ -0.980 and a positive product control (spike) recovery of 50-200%. A coefficient of variation (CV) of less than 25% is expected.

Internally, the PTS reader measures the reaction time in each channel. An archived standard curve specific for each lot of cartridges is constructed using the log of the reaction time vs. the log of the concentration. The sample and spike values are calculated by interpolation of the standard curve using the reaction times. See Certificate of Analysis for the lot specific archived curve range.

TYPES OF ASSAYS

Initial Qualification: Each new lot of cartridges must be qualified upon receipt. The initial qualification testing requires one cartridge with LAL Reagent Water as a sample. The evaluation must demonstrate no detectable endotoxin and acceptable spike recovery (50-200%).

Inhibition/ Enhancement: Absence of interference is demonstrated by achieving acceptable spike recovery (50-200%) on a given sample preparation.

NOTE: A validation currently requires three lots of product to be tested at the same dilution.

Inhibition is usually concentration dependent and can be overcome by dilution with LAL Reagent Water. The most common sources of inhibition are 1) conditions that interfere with the enzyme activity due to ionic strength and/or pH, and 2) those that alter the dispersion of the

endotoxin (positive) control.⁶ If the positive product control fails and a pH related problem is suspected, the pH of the test specimen should be measured to assure a pH within the range of 6.7 - 7.3. Use an endotoxin-free TRIS buffer (Endosafe® BT101, BT103 or equivalent) if pH adjustment is necessary. Do not arbitrarily adjust the pH of unbuffered solutions.

Maximum Valid Dilution: The U.S. Food and Drug Administration has established endotoxin limits of 5 EU/kg for intravenous drugs and 0.2 EU/kg for intrathecal drugs.⁷ Specific limits for compendial items have been adopted.⁹ These limits may be used to determine the extent of dilution that may be used to overcome an interference problem without exceeding the limit endotoxin concentration.¹¹

The Maximum Valid Dilution (MVD) is calculated by formulae presented in the previously mentioned documents and other pharmacopeia.^{7,9}

For drug products that have a published limit, the MVD may be calculated by the following formula:

$$\text{MVD} = \frac{\text{Endotoxin Limit} \times \text{Product Potency}}{\text{Lambda}}$$

EL = K/M, where K= 5EU/Kg and
M= Maximum Dose per kg of
body weight administered per
hour

Product Potency = concentration of
product

Lambda = sensitivity (lowest point on
the archived curve) of test cartridge

For example, the compendial limit for Sterile Water for Irrigation (SWI) is 0.25 EU/mL. If a test cartridge with an archived standard curve containing the lowest level of 0.05 EU/mL of endotoxin is used to test this product, where the potency is 1 mL/mL, the MVD equals 1:5. Thus, SWI may be diluted up to 1:5 to resolve potential inhibition (one part to a total of five parts LRW).

Routine Tests with the Endosafe® - PTS cartridges:

See the User's Guide supplied with the Endosafe® - PTS or the Endosafe® MCS reader for complete operations, procedures, and guidelines.

BIBLIOGRAPHY

1. Bang, F.B., "A Bacterial Disease of *Limulus Polyphemus*", *Bulletin of Johns Hopkins Hospital*, Nr. 98, p. 325 (1956).
2. Cooper, J.F., Levin, J., and Wagner, H.N., "Endotoxins as a Cause of Aseptic Meningitis after Radionuclide Cisternography", *Journal of Nuclear Medicine*, Nr. 16, p. 809 (1976).
3. Cooper, J.F., Levin, J., and Wagner, H.N., "Quantitative Comparison of In Vitro and In Vivo Methods for the Detection of Endotoxin", *Journal of Laboratory and Clinical Medicine*, Nr. 78, p. 138 (1971).
4. Hochstein, H.D., "The LAL Test versus the Rabbit Pyrogen Test for Endotoxin Detection: Update '87", *Pharmaceutical Technology*, Nr. 11(6), p. 124 (1987).
5. Levin, J., and Bang, F.B., "Clottable Protein in *Limulus*: Its Localization and Kinetics of Its Coagulation by Endotoxin", *Thromb. Diath. Haemorrh.*, Nr. 19, p. 186 (1968).
6. Cooper, J.F., "Resolving LAL Test Interferences", *Journal of Parenteral Science and Technology*, Nr. 44:1, p. 13, (1990).
7. Guideline on Validation of the *Limulus Amebocyte Lysate* Test as an End-Product Endotoxin test for Human and Animal Parenteral Drugs, Biological Products, and Medical Devices, U.S. Dept. of Health and Human Services, FDA, December 1987.
8. Interim Guidance for Human and Veterinary Drug Products and Biologicals, U.S. Dept. of Health and Human Services, FDA, July 15, 1991.
9. "Bacterial Endotoxins Test." In The U.S. Pharmacopeia, 26th revision, p. 2023-6, The United States Pharmacopeial Convention, Inc., Rockville, MD., 2002
10. McCullough, K.Z. "Process Control: In-process and Raw Material Testing Using LAL." *Pharmaceutical Technology*, 12(5) p.40 (1988).
11. Weary, M.E. "Understanding and setting endotoxin limits." *Journal of Parenteral Science and Technology*, 44:1, p. 16 (1990).

PATENT INFORMATION

U.S. Patent No: US D472,324 S
U.S. Patent No: US 7,329,538 B2
Other patents pending.

Manufactured By: CHARLES RIVER ENDOSAFE
Div. of Charles River Laboratories, Inc.
1023 Wappoo Road, Suite 43B
Charleston, SC 29407, USA
PHONE NUMBER: 843-766-7575
FAX NUMBER: 843-766-7576

PIPTS20E01